

THE STRENGTH OF TRADITION



# VENOSAN®

2000



Medical Compression Stockings  
and Pantyhose

finest  
*comfort*

# VENOSAN® 2002

## Compression class II

Pressure at the ankle 23–32 mm Hg



Toe	without	without	without	without	without	without	without
	<b>AD</b> Below knee	<b>AF</b> Mid thigh length	<b>AFH</b> Mid thigh length with self supporting plain top	<b>AG</b> Thigh length	<b>AGH</b> Thigh length with self supporting plain top	<b>AGG</b> Thigh length with belt left or right	<b>AT</b> Pantyhose

### long


Article	Colour	Toe	SMALL		MEDIUM		LARGE	
			NORMAL	PLUS	NORMAL	PLUS	NORMAL	PLUS
<b>AD</b>	Bolero	open	SG 54009	SG 54013	SG 54010	SG 54014	SG 54011	SG 54015
<b>AF</b>	Bolero	open	SG 54057	SG 54061	SG 54058	SG 54062	SG 54059	SG 54063
<b>AFH</b>	Bolero	open	SG 54657	SG 54661	SG 54658	SG 54662	SG 54659	SG 54663
<b>AG</b>	Bolero	open	SG 54105	SG 54109	SG 54106	SG 54110	SG 54107	SG 54111
<b>AGH</b>	Bolero	open	SG 54605	SG 54609	SG 54606	SG 54610	SG 54607	SG 54611
<b>AGG</b>	left	open	SG 54153	SG 54157	SG 54154	SG 54158	SG 54155	SG 54159
	right	open	SG 54185	SG 54189	SG 54186	SG 54190	SG 54187	SG 54191
<b>AT</b>	Bolero	open	SG 54217	SG 54221	SG 54218	SG 54222	SG 54219	SG 54223

### short

Article	Colour	Toe	SMALL		MEDIUM		LARGE	
			NORMAL	PLUS	NORMAL	PLUS	NORMAL	PLUS
<b>AD</b>	Bolero	offen	SG 54001	SG 54005	SG 54002	SG 54006	SG 54003	SG 54007
<b>AF</b>	Bolero	offen	SG 54049	SG 54053	SG 54050	SG 54054	SG 54051	SG 54055
<b>AFH</b>	Bolero	offen	SG 54649	SG 54653	SG 54650	SG 54654	SG 54651	SG 54655
<b>AGG</b>	left	offen	SG 54145	SG 54149	SG 54146	SG 54150	SG 54147	SG 54151
	right	offen	SG 54177	SG 54181	SG 54178	SG 54182	SG 54179	SG 54183
<b>AT</b>	Bolero	offen	SG 54209	SG 54213	SG 54210	SG 54214	SG 54211	SG 54215

# VENOSAN® 2000

## Seamless made-to-measure stockings

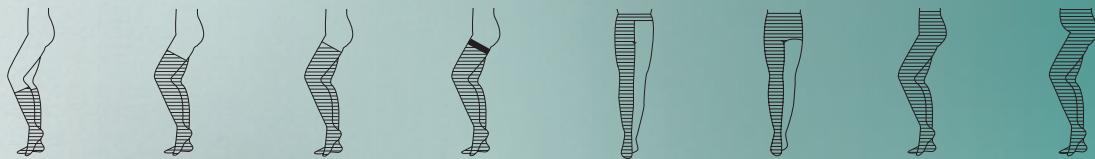

Toe	with /without	with /without	with /without	with /without	with /without	with /without	with /without	with /without
	<b>AD</b> Below knee	<b>AF</b> Mid thigh length	<b>AG</b> Thigh length	<b>AGH</b> Thigh length with self supporting top	<b>AGG</b> Thigh length with belt left or right	<b>AGT</b> One-leg-Panty	<b>AT</b> Pantyhose	<b>ATU Materna</b> Maternity Pantyhose

**2002**


Ccl. II

**2003**

Ccl. III

**2004**

Ccl. IV

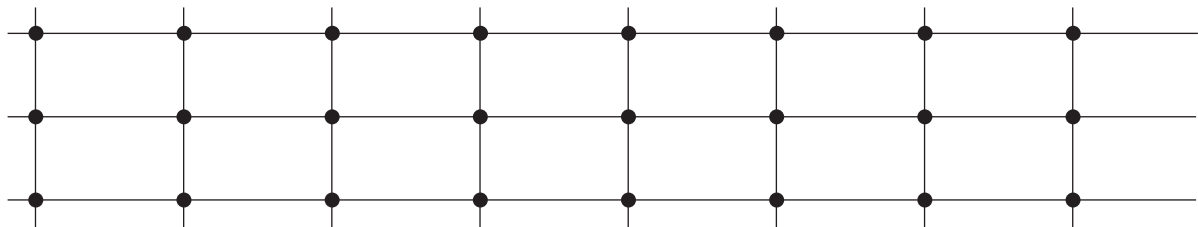

Measurements should be taken with the patient standing and if possible, after the leg has been rested and swelling is at a minimum.

All measurements should be made on the bare leg.

Custom made stockings may not be returned for credit.

The manufacturer cannot be held responsible for incorrect measurements.

**Please ask for the special order form for made-to-measure stockings.**

# VENOSAN® 2003

Compression class III

Pressure at the ankle 34–46 mm Hg



Toe	without	without	without	without	without	without	without
	<b>AD</b> Below knee	<b>AF</b> Mid thigh length	<b>AFH</b> Mid thigh length with self supporting plain top	<b>AG</b> Thigh length	<b>AGH</b> Thigh length with self supporting plain top	<b>AGG</b> Thigh length with belt left or right	<b>AT</b> Pantyhose

## long

Article	Colour	Toe	SMALL		MEDIUM		LARGE	
			NORMAL	PLUS	NORMAL	PLUS	NORMAL	PLUS
<b>AD</b>	Bolero	offen	SG 55009	SG 55013	SG 55010	SG 55014	SG 55011	SG 55015
<b>AF</b>	Bolero	offen	SG 55057	SG 55061	SG 55058	SG 55062	SG 55059	SG 55063
<b>AFH</b>	Bolero	offen	SG 55657	SG 55661	SG 55658	SG 55662	SG 55659	SG 55663
<b>AG</b>	Bolero	offen	SG 55105	SG 55109	SG 55106	SG 55110	SG 55107	SG 55111
<b>AGH</b>	Bolero	offen	SG 55605	SG 55609	SG 55606	SG 55610	SG 55607	SG 55611
<b>AGG</b>	left	offen	SG 55153	SG 55157	SG 55154	SG 55158	SG 55155	SG 55159
	right	offen	SG 55185	SG 55189	SG 55186	SG 55190	SG 55187	SG 55191
<b>AT</b>	Bolero	offen	SG 55217	SG 55221	SG 55218	SG 55222	SG 55219	SG 55223

## short

Article	Colour	Toe	SMALL		MEDIUM		LARGE	
			NORMAL	PLUS	NORMAL	PLUS	NORMAL	PLUS
<b>AD</b>	Bolero	offen	SG 55001	SG 55005	SG 55002	SG 55006	SG 55003	SG 55007
<b>AF</b>	Bolero	offen	SG 55049	SG 55053	SG 55050	SG 55054	SG 55051	SG 55055
<b>AFH</b>	Bolero	offen	SG 55649	SG 55653	SG 55650	SG 55654	SG 55651	SG 55655
<b>AGG</b>	left	offen	SG 55145	SG 55149	SG 55146	SG 55150	SG 55147	SG 55151
	right	offen	SG 55177	SG 55181	SG 55178	SG 55182	SG 55179	SG 55183
<b>AT</b>	Bolero	offen	SG 55209	SG 55213	SG 55210	SG 55214	SG 55211	SG 55215



### Material

60% Nylon Polyamide

25% Lycra® Elastane

15% Cotton

100% Cotton next to the skin

Bolero

# VENOSAN®

How should you take measurements?

For correct measurements it is very important to take your time and proceed systematically. Wherever possible, the measurements should be taken when the leg is rested and not swollen. **The measurements should in any case be taken with the patient standing.** Well fitting medical compression stockings are the key to a successful treatment and satisfied patients.



AD



AF



AG



AGG



AT

Procedure:

Measure:

### Below knee AD

1. the length of the leg from the base of the heel to just below the knee **ID**
2. the circumference of the ankle **cB**
3. the circumference of the calf **cC**
4. the circumference of the leg just below the knee **cD**

Measure:


### Mid thigh length AF, AFH

1. the length of the leg from the base of the heel to the middle of the thigh **IF**
2. the circumference of the ankle **cB**
3. the circumference of the calf **cC**
4. the circumference of the leg just below the knee **cD**
5. the circumference of the leg in the middle of the thigh **cF**

Measure:

### Thigh length AG, AGH, thigh length with belt AGG and Pantyhose AT

1. the inside leg **IG** (for AGG/AT up to the groin **IK**)
2. the circumference of the ankle **cB**
3. the circumference of the calf **cC**
4. the circumference of the leg just below the knee **cD**
5. the circumference of the leg in the middle of the thigh **cF**
6. the circumference of the thigh at the widest point, usually about 5 cm below the groin **cG**

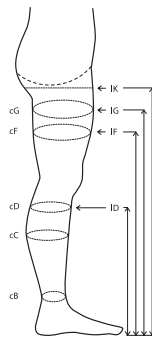

## MEASUREMENT POSITIONS / CIRCUMFERENCE MEASUREMENTS

**cG** Thigh (5 cm below groin)

**cF** Middle of the thigh

**cD** Narrowest point below the knee

**cC** Widest part of the calf

**cB** Ankle (narrowest point above ankle bone)

Subsequently refer to the size chart below to see:

1. the length of the leg **ID**, **IF**, **IG** or **IK** (short or long)
2. in which size column the measured circumference of the ankle **cB** fits (small, medium or large)
3. in which size column the other circumferences **cC**, **cD**, **cF** and **cG** fit (normal or plus)

You have now selected the correct size, eg. medium plus short and you will find the corresponding article number on the article chart.

## SIZE CHART

### AD

Circumference (cm)	SMALL		MEDIUM		LARGE	
	NORMAL	PLUS	NORMAL	PLUS	NORMAL	PLUS
<b>cB</b>	19 – 23	19 – 23	23 – 26	23 – 26	26 – 29	26 – 29
<b>cC</b>	30 – 36	35 – 40	35 – 40	39 – 45	39 – 45	44 – 49
<b>cD</b>	28 – 34	33 – 38	33 – 38	37 – 43	37 – 43	42 – 47

NORMAL size = for standard legs

PLUS size = for patients with a wide calf and large below knee circumference

### AF/AG/AGG/AT

Circumference (cm)	SMALL		MEDIUM		LARGE	
	NORMAL	PLUS	NORMAL	PLUS	NORMAL	PLUS
<b>cB</b>	19 – 23	19 – 23	23 – 26	23 – 26	26 – 29	26 – 29
<b>cC</b>	30 – 36	30 – 36	35 – 40	35 – 40	39 – 45	39 – 45
<b>cD</b>	28 – 34	28 – 34	33 – 38	33 – 38	37 – 43	37 – 43
<b>cF</b>	42 – 50	48 – 56	48 – 56	54 – 62	54 – 62	60 – 71
<b>cG</b>	44 – 52	53 – 62	50 – 58	60 – 69	56 – 66	67 – 77

NORMAL size = for standard legs

PLUS size = for patients with wide thighs

Length measurement (cm)	KURZ	LANG
<b>ID</b> AD	36 – 40	40 – 44
<b>IF</b> AF	56 – 62	62 – 68
<b>IG</b> AG		72 – 79
<b>IK</b> AGG/AT	67 – 74	74 – 85

*Salzmann*  
**MEDICO**

SALZMANN AG

Salzmann MEDICO

Rorschacher Strasse 304

CH-9016 St.Gallen

Tel. (+41) 071 282 12 12

Fax (+41) 071 282 12 10

medico.sg@salzmann-group.ch

www.salzmann-group.ch

www.venosan.com

